

# i3 ONE sterile filter







# **INTENDED USE & FIELDS OF APPLICATION**

The i3 ONE point-of-use sterile filter is used in medical settings to prevent infections caused by bacteria such as *Legionella pneumophila*, *Pseudomonas aeruginosa* and mycobacteria. It is used on the water faucet for the sterile filtration of tap water. The water filter protects patients receiving medical treatment, for example in oncological, transplant and intensive care wards, and also neonates, from infections caused by waterborne pathogens.

# **FEATURES**

CE-marked Class I medical device

Bacterial retention as per ASTM F838

Sterile filtration with a nominal pore size of 0.2 µm (microfiltration)

Duration of use: 50 days (7 weeks)

High, stable flow rate throughout the period of use

Protection due to the bacteriostatic properties of the filter housing in accordance with ISO 22196

Proven biocompatibility to protect the user

Filter surface can be disinfected by wiping

Compact design means use of the water tap is not limited

Angled laminar water jet

Connected to water tap by means of quick-release coupling

Clear traceability

# **TECHNICAL SPECIFICATIONS**



## **APPLICATIONS**

Water tap



## **FLOW RATE**

12 l/min at 5 bar



# **AUTHORIZED FOR USE** WITH DRINKING WATER

Plastics suitable for use with drinking

water as per BS 6920



Medical



## **BACTERIAL RETENTION**

> 7 log (Brevundimonas diminuta) per cm<sup>2</sup> of filter area



#### **DISINFECTION**

in 24 hr.

Can be disinfected by wiping with normal surface disinfectant

**RETROGRADE PROTECTION** 

> 4 log Pseudomonas aeruginosa,

Escherichia coli, Staphylococcus aureus



#### **OPERATING TEMPERATURE**

60 °C (max 70 °C for 30 min)



## **DURATION OF USE**

50 days (7 weeks)



#### **OPERATING PRESSURE**

max 5 bar



## **SHELF LIFE**

5 years



#### **MEMBRANE**

i3 high performance flat membrane



## **MATERIAL AUTHORIZATION**

Product is biocompatible as per ISO 10993-1



## **CHLORINE RESISTANCE**

≤ 200 ppm active chlorine for 1 hour ≤ 10 ppm active chlorine for 50 days



## **ATTACHMENT**


Quick-release coupling is required for installation



# **LABELLING**

Serial number / DataMatrix Code for documentation of filter changes

# **DIMENSIONS (mm)**

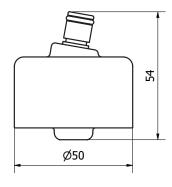

# **ORDER INFORMATION & ACCESSORIES**



**i3 ONE** Sterile filter, water tap 203100 (Item no.)



### Packaging units 2 pcs / 12 pcs



M24 guick-release coupling (OT) 2060001 (Item no.)



**M22** quick-release coupling (IT) 2060000 (Item no.)









